

HIPRO

HIGH INTENSITY FOCUSED ULTRASONIC



HIPRO

HIGH INTENSITY FOCUSED ULTRASONIC

What is HIPRO?

HIPRO delivers HIFU energy to SMAS layer (which is located between subcutaneous tissues and muscle) and makes TCP(Thermal Coagulation Points) zone without damaging skin surface by 4.5mm cartridge, also by using 3.0mm and 1.5mm cartridges, the HIFU energy can stimulate the collagen to be remodeled for skin tightening.

Characteristics of HIPRO



- CE & KFDA approval as total face lifting, tightening & body slimming
- Five types of cartridge 1.5mm, 3.0mm, 4.5mm, 8.0mm and 13.0mm
- The shortest treatment time by continuous mode
- Reasonable price of consumables
- Less pain than other HIFU device



Treatments of HIPRO for lifting & tightening




- Forehead wrinkle lifting
- Nasolabial fold
- Wrinkle around the eyes
- Neck Wrinkle




- Abdomen
- Arm
- Thigh
- Back

History of Skin Lifting

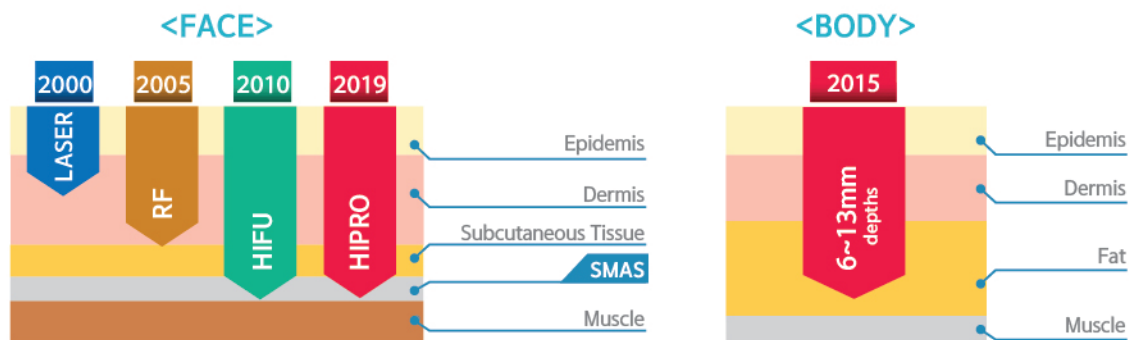

HIFU system is used to make the thermal coagulation to SMAS layer or deep dermis as non-invasive skin lifting and tightening treatment. Now, HIFU system is developed to area of Body for slimming. HIPRO can make thermal coagulation to layer of fat and help to reduce fat by using principle of HIFU Energy.

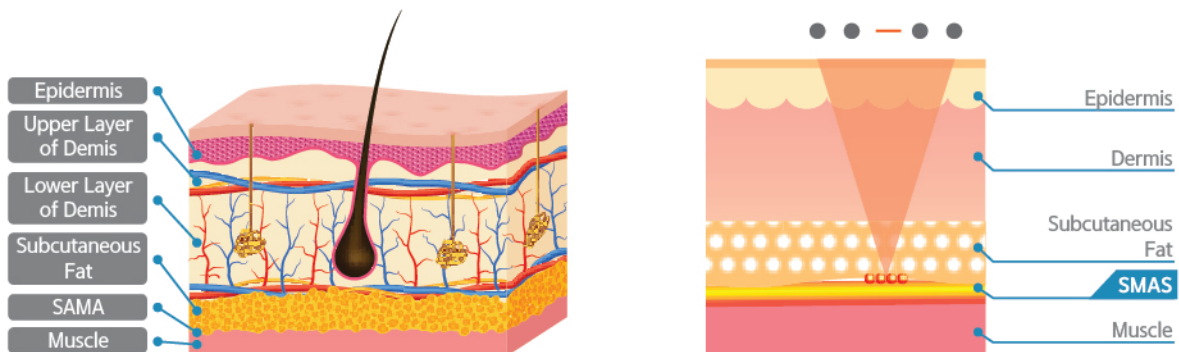

- Make over 60°C of TCP zone on SMAS layer by delivering HIFU energy.
- By wound healing process on TCP zone of SMAS, the surface of muscle contraction, regrowth of collagen and elastin are stimulated by 3.0mm & 1.5mm cartridges.
- Skin lifting and tightening occurs by collagen regeneration process of dermis and the contraction of SMAS layer.

SPECIFICATIONS

Energy Type	HIFU (High Intensity Focused Ultrasonic)
Energy Output	0.1 ~ 2.0 [J]
Input Power	AC 100-230 [V], 50/60 [Hz]
Interval	1.0 ~ 2.0 [mm]
Shot length	30mm(full length), 15mm(half length), dual(full length)
Cartridge A (1.5mm)	10,000 shots
Cartridge B (3.0mm)	10,000 shots
Cartridge C (4.5mm)	10,000 shots
Cartridge D (8mm)	10,000 shots
Cartridge E (13mm)	10,000 shots
Dimension	430 X 570 X 1050 [mm](W x D x H) 430 X 410 X 246 [mm](W x D x H)
Weight	HIPRO : 35 [kg] HIPRO-S : 8 [kg]



Tel : +82.2.2208.0905 E-mail : globalsales@daejumedi.co.kr

www.infilux.co.kr